




CELL-DYN EMERALD 22*

So small ... so great optical

UNI-Flow technology with patented CELL-DYN Emerald 22*
lyse to provide a true optical 5-part WBC differential

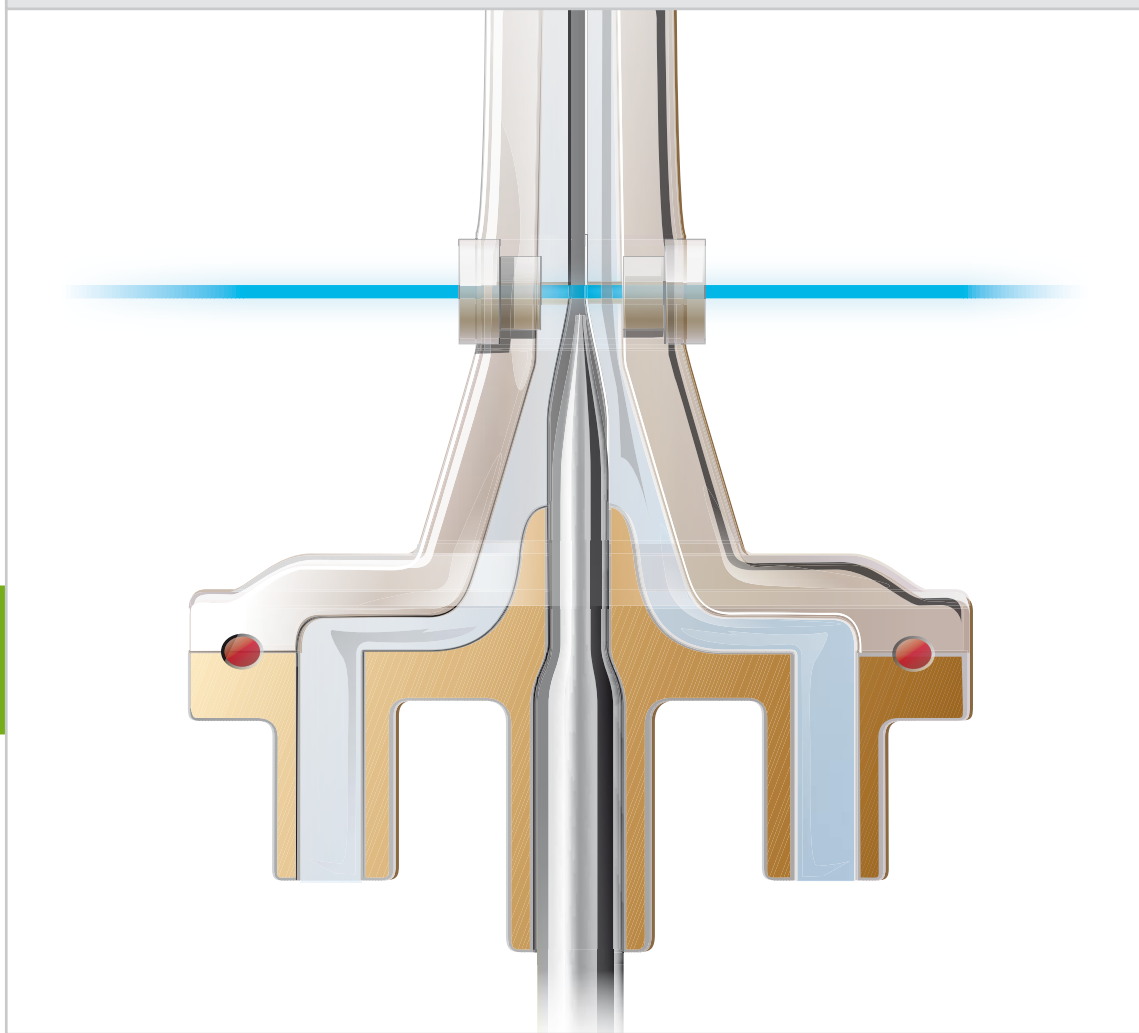

* in development

CELL-DYN
Emerald 22*



Put science on your side.


Abbott
A Promise for Life

A true optical 5 part WBC differential based on UNI-Flow technology

CELL-DYN Emerald 22* uses UNI-Flow technology, which includes a patented cyanide-free lyse reagent and optical bench with flow-cell and blue LED at 455 nm. The unique flow cell design utilizes a passive sheath and a sample stream under pressure, which allows a more diluted sample to be used for analysis.

The CELL-DYN Emerald 22* lyse stabilizes the WBC populations, is used to count the white blood cells (WBC) and measure hemoglobin (Hb) in the same counting chamber at a single dilution.

CELL-DYN Emerald 22* Lyse is ready to use, cyanide-free, and formaldehyde-free!

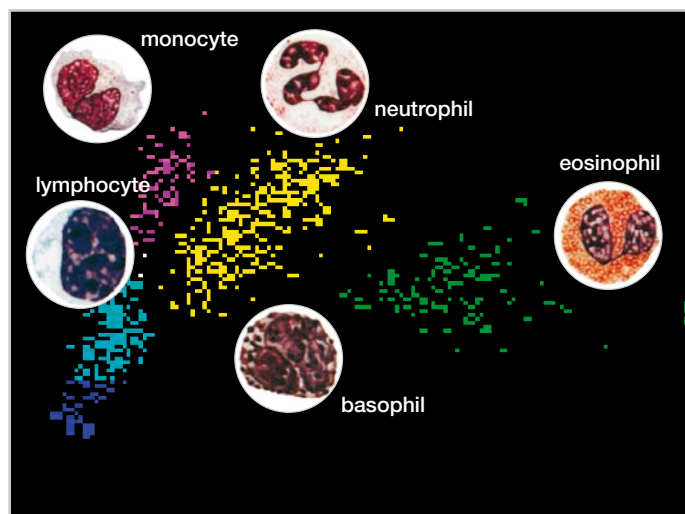

UNI-Flow technology ... so optical

Technology	DC Electrical Impedance (for WBC, RBC, MCV, PLT and MPV) Optical scatter (for WBC differential) Absorption photometry		
Throughput	>45 samples per hour CBC Diff and CBC mode		
Sample Size	17 µl		
Parameters	White Cells WBC Lym # Lym % Neut # Neut % Eo # Eo % Baso # Baso % Mono # Mono %	Red Cells RBC MCH HGB MCHC HCT RDW MCV	Platelets PLT MPV
Physical Dimensions	Height 13.8" (35 cm) Width 9.8" (25 cm) Depth 13.8" (35 cm) Weight –26.5 lbs (12 kg) (without on-board reagents)		